
ORBIS

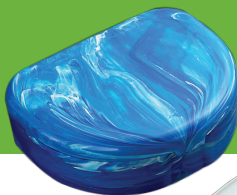StressGuard

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ORBIS

StressGuard

The comfortable, easy-fit
night-time dental guard

Protects teeth from the effects of grinding and bruxism

Designed by an Orthodontist

Orbis StressGuard night-time dental guard has a comfortable, easy-fit design to protect your teeth from the effects of night-time grinding (bruxism).

Not recommended for patients under 18 years of age

Not suitable as athletic mouth guard.

Not recommended for use on patients with:

- ♪ loose teeth
- ♪ TMJ (temporo mandibular joint disorder)
- ♪ any tooth or jaw pain

Checklist for patient:

- ♪ no loose fillings or crowns
- ♪ cavities with no fillings
- ♪ no pain in jaws, teeth or face or difficulty in chewing
- ♪ no missing teeth or mouth sores
- ♪ no gum disease
- ♪ no serious breathing, respiratory or other health problems

When using this product you are advised to visit your dentist every six months

Instruction for use

Place the Orbis StressGuard between the upper teeth and lip and the bite tabs between your upper and lower teeth.

Can also be used on lower teeth.

Orbis StressGuard is designed to keep your lower and upper teeth apart

It does not fit tightly around your teeth.

It can be trimmed with scissors or silicone trimmer.

We recommend that you rinse it with soap and water after use.

When first using the Orbis StressGuard:

- ♪ you are advised to wear it for a couple of hours to get used to it
- ♪ more salivation may occur which is normal