



MIDAS



CHAIRSIDE RESTORATIONS HAVE NEVER BEEN

SMART **AND** SIMPLE

Stress-free, chairside.

The world's first digital restoration press and the simplest and smartest method of 3D printing dental restorations.

Patented Digital Press Stereolithography (DPS)

Midas is an all-new solution for chairside restorations. Digital Press Stereolithography combines the tank, resin, and platform into a single capsule, allowing it to print with materials that are more viscous than conventional 3D printers.

Print a range of restorations in 10 minutes



Up to **3** OR Up to **6** OR Up to **9**
Crowns **Inlays/Onlays** **Veneers**

Combination of applications across all 3 capsules

Midas Capsule

Midas accommodates up to three capsules at once





Validated Materials

Features SprintRay Ceramic Crown, Class II FDA and MDR cleared for use in 3D printing definitive individual and fixed full single crowns, definitive partial crowns in the front and posterior area, individual and fixed single veneers. Prints single-unit crowns, veneers, inlays, or onlays in less than 10 minutes.

Effortless Fabrication

Just scan, place, and press – it's as simple as that. Experience no-mess handling with fully enclosed resin capsules. Ready to use right out of the box, with no extensive training required.



Step 1 Scan

Scan the capsule using the barcode scanner.



Step 2 Place

Place the capsule where the computer instructs.



Step 3 Press

Midas will press up to 3 capsules in 10 minutes.



Step 4 Remove

Remove the build platform – your restoration is ready.


Design In Minutes With AI Studio

AI Studio combines the convenience of AI design with the versatility of CAD to create a simple, cloud-based design suite. Designing a restoration with Studio is easy for anyone, while customization options like margin and adaptation adjustments help create a perfect fit.


When you buy Midas, you'll unlock access to AI Studio and unlimited AI designs for one year.




Upload Scans

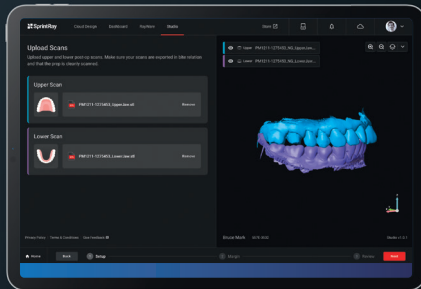

Confirm Margins

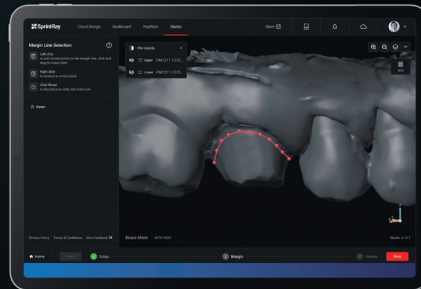

Final Adjustments

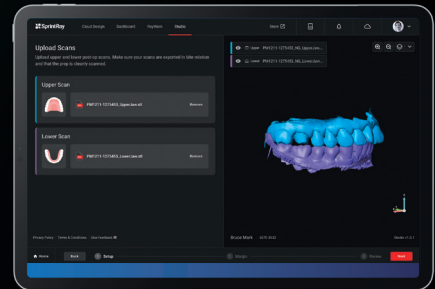

Midas Restorative Workflows are supported by our trusted Dental Partners:

 **DEXIS** **exocad** **iTero™** **ivoclar** **MEDIT** **3shape** 



SprintRay Europe GmbH
Brunnenweg 11
64331 Weiterstadt | Germany
www.sprintray.com
info.eu@sprintray.com



See Midas in Action.
Request a demo today!