

Give the gift of a beaming smile





Proxeo: the dental prophylaxis for a beaming smile

Professional tooth cleaning is an effective method of preventing plaque deposits from forming, giving your patients a beaming smile. With Proxeo, W&H offers a compact, easy to handle solution, which makes it even easier to give your patients what they want.

Proxeo has two systems for everything you need, no matter whether it is cleaning, polishing or fluoridating!

- > UNIVERSAL system

 Easily combined with all common cups and brushes
- > YOUNG system Can be used with all YOUNG cups and brushes



The advantages are clear:

- > Flexible and fast
- > Long lifespan as a result of the unique sealing systems
- > Low-maintenance and hygiene-friendly
- > The small head gives full visibility of the area to be treated
- > Efficient cleaning assured by using optimum speeds
- > Ergonomic shape

Exceptionally adaptable: the Universal system

Change in the blink of an eye

Very easy to handle: Thanks to the screw-in or snap-on technology, you are free to screw on or attach all your cups and brushes however you want.

Longer satisfaction

With this sealing system specially developed by W&H, paste and treatment particles stay where they belong and do not penetrate into the instrument. This guarantees a long lifespan and total satisfaction with every use.






Versatile and flexible: the YOUNG system

The YOUNG Hygiene Plus

The YOUNG disposable contra-angle handpiece can easily be fitted onto the Proxeo handpiece. The entire contra-angle handpiece is removed after use. It is extremely user-friendly and facilitates maintenance.

On top form

The YOUNG system's unique Triple Seal sealing prevents even the smallest particles of dirt from entering the instrument. This seal is renewed every time the cup is changed. This ensures the highest possible level of hygiene and safety. It also helps preserve the instruments and therefore ensures longevity.

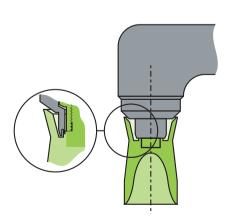



Whether it's YOUNG or UNIVERSAL: Proxeo outshines everything else

All clean

Proxeo is thermo washer disinfectable and sterilizable up to 135°C.





Superior quality guaranteed

In addition to their precision, the progressive technology and the unrivalled quality, all Proxeo instruments also come with a 24-month warranty - which means that you can expect the very best from our products in every respect.

Improved visibility

There is nothing in the way: the small contraangle head means that you will always have everything in view - even treatment sites in the distal area.

Optimum speed

The rotational speed and the transition ratio of 4:1 enable:

- > greater cleaning efficiency
- > gentle polishing of tooth surfaces
- > no splattering of paste



Cleaning and polishing of tooth surfaces in the interdental region (e. g. using the snapon system)



Polishing of composite fillings in the frontal region (e. g. using screw-in system)



Polishing of amalgam fillings in the interdental region (e. g. using screw-in system)

Technical data

C €	Universal system Once of the state of the s	Young sys	PICKED
Туре:	WP-64 MU	WP-64 M	HP-44 M
Max. motor speed:	10,000 rpm	10,000 rpm	10,000 rpm
Transition ratio:	4:1	4:1	4:1
Recommended speed for cups and brushes:	1,250 1,500 rpm	1,250 1,500 rpm	1,250 1,500 rpm
Maximum speed for cups and brushes:	2,500 rpm	2,500 rpm	2,500 rpm
Attachments that can be used:	All common screw-in cups and brushes [1-72 UNF 2A screw thread] as well as snap-on cups and brushes [adaptor enclosed]	YOUNG screw-in cups and brushes	Disposable contra-angle handpiece with plastic shaft, Doriot system
Coupling system:	For motors with ISO connection in accordance with ISO 3964/DIN 13940	For motors with ISO connection in accordance with ISO 3964/DIN 13940	For motors with ISO connection in accordance with ISO 3964/DIN 13940
Length:	89 mm	89 mm	64 mm
Weight:	71 g	69 g	57 g
Sterilizable:	In the autoclave up to 135°C	In the autoclave up to 135°C	In the autoclave up to 135°C

W&H Dentalwerk Bürmoos GmbH Ignaz-Glaser-Straße 53, Postfach 1 5111 Bürmoos, Austria t +43 (0) 6274 6236-0 f +43 (0) 6274 6236-55 office@wh.com