

U600

Ultrasonic Scaler

Pleasant & Efficient




Intelligent power output
Experience painless treatment

U600 Ultrasonic Scaler

A global professional scaler manufacturer focuses on technological innovation for 30 years.

International high-tech enterprise
30 years innovation
Woodpecker scaler owns **101** national patents
9 national invention patents

Effective and thorough endodontic irrigation

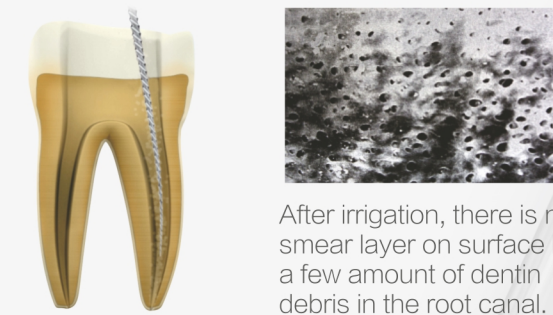

After irrigation, there is no smear layer on surface but a few amount of dentin debris in the root canal.


Root canal irrigation is the most widely used clinical support methods. Its purpose is to remove the substance inside root canal, dissolve diseased tissue, remove pathogenic microorganisms, and remove smear layer. Lubricate the root canal internal wall. During root canal preparation process, prevent the infectious substance from being pushed to the deep part or out of the apical foramen, reducing the risk of instrument fracture in root canal.

Compared with the traditional syringe irrigation, ultrasonic irrigation can better remove bacteria, dentin debris, and smear layer in root canal. And it can well clean special structures such as the lateral branch canal and root canal narrow area. Currently, ultrasonic irrigation is widely used in endodontic therapy in order to receive a better root canal preparation effect.

Sine wave drive

Avoid instantaneous stimulation caused by power surge, ensuring comfortable treatment.

U600 Mild power output

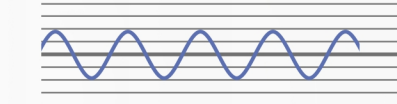


Rectangular power out of other brand scaler




High-quality handpiece

Aluminum shell can be autoclaved, not easy to cause cracking or deformation. Latest handpiece sealing technique enables long durability. User-friendly design brings better hand feeling.

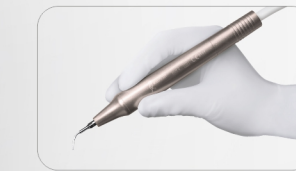

Automatic water supply system

Easy to change water, eliminating redundant steps and increasing the working efficiency. Solutions such as hydrogen peroxide, sodium hypochlorite, and chlorhexidine can be used. Clear 50ml scale mark makes it convenient to formulate solutions.

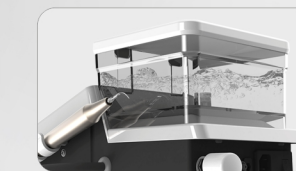

> High configuration

G4 x2	G12 x2	P1 x2	P3 x2	P8 x2	E60 x2	E62 x2
----------	-----------	----------	----------	----------	-----------	-----------

Standard configuration


G4 x1	G12 x1	P1 x1	P3 x1	P8 x1	E60 x1	E62 x1
----------	-----------	----------	----------	----------	-----------	-----------

Note: Woodpecker reserves the right of modifying the tip configuration without further notice.

Supragingival scaling





Implant maintenance

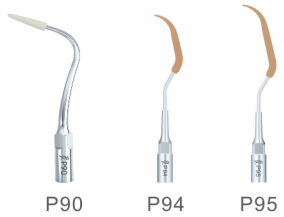

Subgingival scaling



Root canal treatment

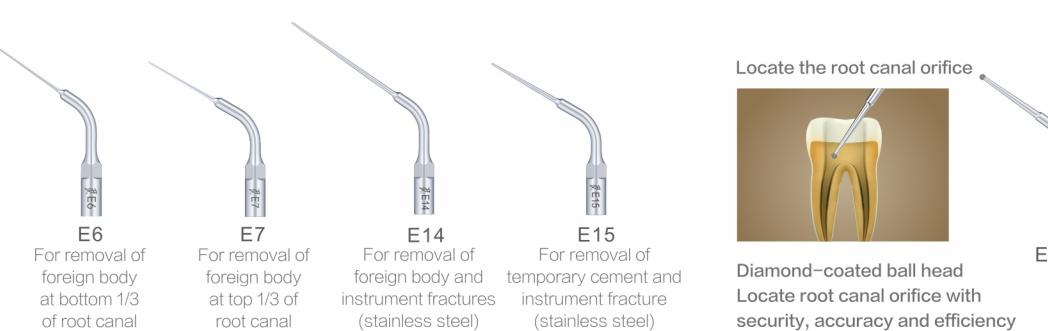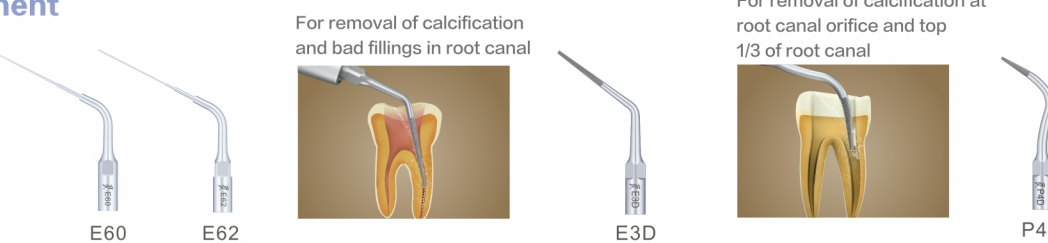

Cavity preparation



GUILIN WOODPECKER MEDICAL INSTRUMENT CO., LTD.

Registration Add & Manufactory Add:
Information Industrial Park,Guilin National High-Tech
Information Park 541004 Guilin,Guangxi China
Tel:
Europe Sales Dept.: +86-773-5873196, +86-773-2125222
North America, South America &
Oceania Sales Dept.: +86-773-5873198, +86-773-2125123
Asia & Africa Sales Dept.: +86-773-5855350, +86-773-2125896
Fax: +86-773-5822450
E-mail: woodpecker@glwoodpecker.com, sales@glwoodpecker.com
Website: http://www.glwoodpecker.com
ZMN/WI-10-099 V1.0

Distributor:



Scan and login website
for more information

