


BROCHURE

DEKO Modular Rack

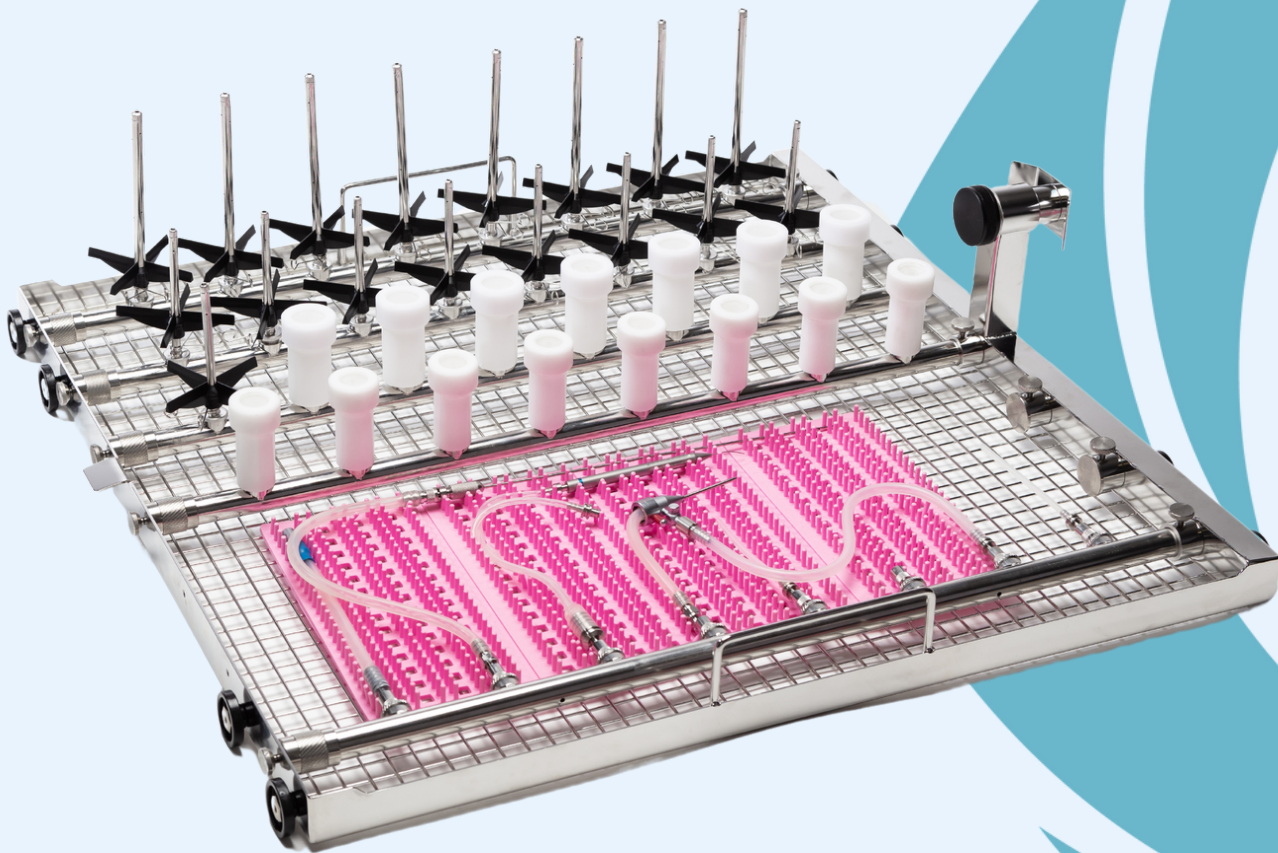

Preventing infections since 1976.

DEKO
MEDTECH

Modular Rack

Experience streamlined workflows and optimal efficiency with the **DEKO Modular Rack concept**. Whether processing dental or MIS instruments, glassware or something totally different, the Modular Rack easily adapts to specific needs.

Tailored solutions for unique needs

Available for **DEKO 260, DEKO 2000 and DEKO D32**, the Modular Rack offers custom solutions for the most challenging, versatile and unique needs our customers have. A wide array of washing accessories, adapters and inserts for fully customizable washing configurations.

Whether the specific washing needs are limited or require extensive capacity, with the Modular Rack it is easy to form the perfect solution for each occasion.

Limited space, mixed loads? The Modular Rack is the answer

Space is often a limited resource. With the adaptable DEKO Modular Rack **efficiency and quality go hand-in-hand** and there is no need for compromises.

The rack provides convenience with easy adaptability to various washing needs with simple toolless solutions replacing bulky specialty racks taking up valuable space.

Mixed loads and frequent changes in load compositions? Not a problem for the Modular Rack! Select the washing accessories fit for your specific needs and adjust the rack in just a few minutes.



Capacity

The Modular Racks provide impressive capacity. Equip the rack base with the desired number of branches and washing accessories that best fit your needs. Up to 90 (ninety) instrument couplings 100% adaptable with a range of accessories.

DEKO 260	6 branches	42 couplings
DEKO 2000	7 branches	56 couplings
DEKO D32	9 branches	90 couplings

Capacity
that makes
a difference



Adaptability made simple

Solutions that work


All the washing accessories can be easily replaced without tools. Choose from a selection of spray jets, instrument holders, luer locks and various other accessories to **curate the perfect combination** to suit your washing needs.

Rack branches can be secured at three **different angles** to accommodate various instruments.

Water is introduced to the instruments through **a fine mesh filter**. DEKO D32C (B) Modular Rack base is equipped with double inlets for maximized washing efficiency.

Adding trays, basket or inserts on the modular base allow for practically **unlimited rack designs**.

Modular Rack concept

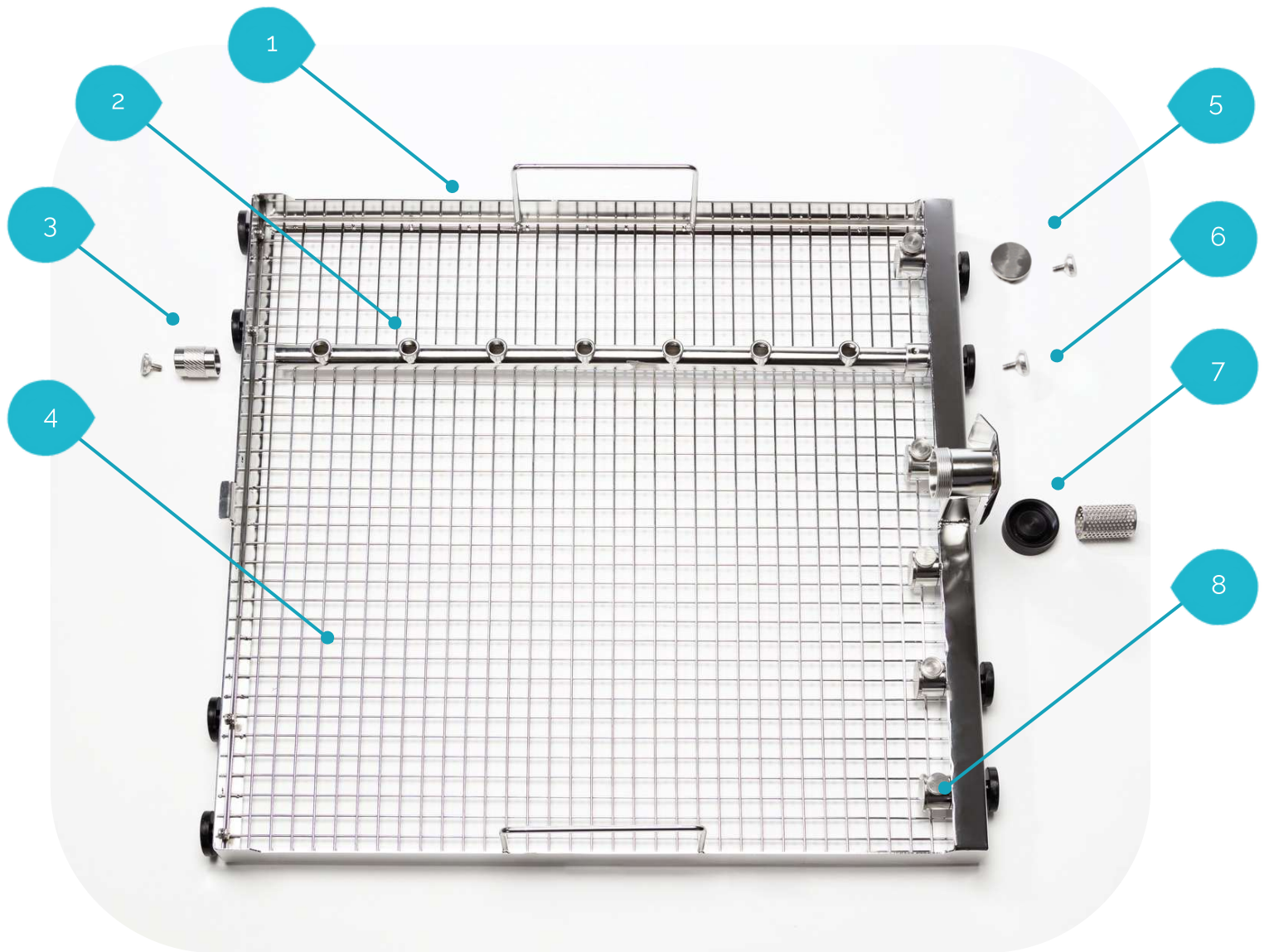

1

Modular Rack base

2

Rack branch with accessories couplings

3

Plug for branch tube + fixing screw for plug

4

Space for trays & inserts

Plug for muff / branch coupling + fixing screw for plug

5

Fixing screw to fasten branch in desired angle

6

Fine mesh filter + filter cap

7

Plugged branch coupling

8


Modular Rack accessories - right fit for every need

BRANCHES

Branch for DEKO 260

Branch for DEKO 2000

Branch for DEKO D32

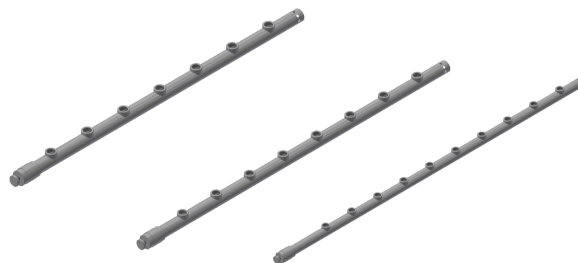

SPRAY JETS

Spray jet Ø3x100

Spray jets Ø4x60, 4x95, 4x125, 4x200

Spray jets Ø6x15, 6x55, 6x100, 6x150, 6x200

Spray jets Ø8x150, 8x200, 8x295, 8x350,

8x350 for respiratory devices




SPRAY JET ACCESSORIES

Bottle rests for Ø4/6/8 jets

Height adjustment springs for Ø4/6/8 jets

Bottle supports for Ø4/6/8 jets

Retainer clips for Ø4/6/8 jets

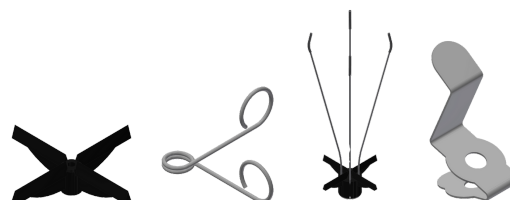

LUER LOCKS

Luer locks w/M12 thread

Luer locks attachable to silicone hose




INSTRUMENT HOLDERS

Instrument holders G3/4", G7/8"

Funnel for dental instruments w/spring

Funnels for pipettes

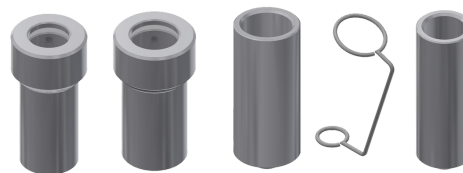


MISC ACCESSORIES

Silicone hoses Ø4/7, Ø17/24, Ø22/28

Nozzle tubes for hoses

Plugs and screws

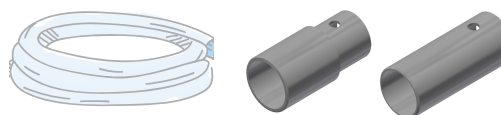

Reference solutions

DEKO 2000 Modular Rack for MIS instruments

- Modular Rack frame for DEKO 2000
- 5 x Rack branches DEKO 2000
- 2 x Plugs+screws for branch connections
- 8 x Spray jets Ø8x200
- 9 x Spray jets Ø6x150
- Bottle rests + height adjustments springs
- 7 x Instrument holders G7/8"
- 8 x Instrument holders G3/4"
- Male luer locks
- Female luer locks
- Silicone hose Ø4/7

Additional items are photo props.



DEKO 260 Modular Rack for dental instruments

- Modular Rack frame for DEKO 260
- 3 x Rack branches for DEKO 260
- 3 x Plugs+screws for branch connections
- 9 x Spray jets Ø6X150
- Bottle rests + height adjustment springs
- 7 x Instrument holders G7/8"
- 7 x Instrument holders G3/4"
- Tray support for 8 trays
- Small fine mesh basket with lid

Additional items are photo props.



Hints and tips



Modular rack base and all the accessories are available separately enabling unique, environmentally conscious and cost-efficient solutions.



The rack base includes plugs and screws for plugging the brach couplings. Rack branches are sold with required plugs & screws for fastening.



Switch effortlessly between purposes with additional branches equipped with pre-installed washing accessories for quick and easy adaptation.



Diversify the washing options with trays and inserts, available on DEKO Racks and Accessories catalogue.



See DEKO Racks and Accessories catalogue for compatibility of washing accessories between standard DEKO racks and the Modular racks.



With separately sold accessories, the customizable rack concept provides a flexible and scalable solution that adapts to changing requirements ensuring long-time value.

Notes



Preventing infections since 1976.

DEKO MedTech Oy | Vartiokuja 1 | FI-76850 Naarajärvi | Finland
www.dekomedtech.com