




Rod attachments

Vario-Soft 3

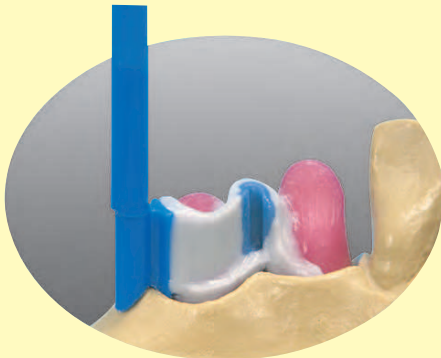


Soft matrices

Matrixes that have proved their reliability for 20 years provide safety and ensure high comfort of wear for the patient.

Soft soft matrices

Special soft plastic compensates small divergences and minor processing imperfections.

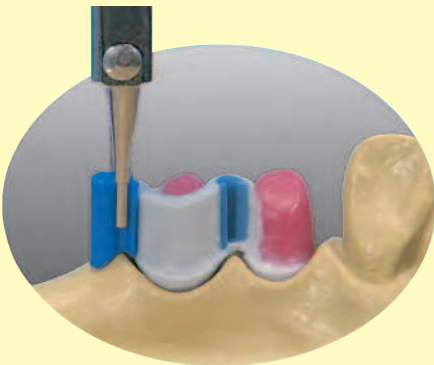


After casting, patrices must only be processed using rubber polishers and high-lustre buffs.




The white duplicating matrix that has been adapted from the basal direction provides the perfect precondition for all other types of friction.



Master model prepared for duplicating.

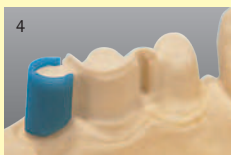

Wax matrix housing on the investment material model guarantees a uniform chrome cobalt housing.

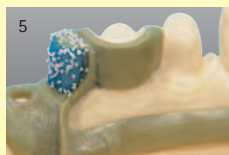

Completed wax pattern of the later chrome cobalt framework.



The use of the inserting instrument ensures precise positioning of the matrices.

Assortment

13 pieces

Vario-Soft 3

2 vs 3 Patrices

1 Matrix inserting instrument

2 Duplicating matrix

2 Wax matrix housing

2 Soft Matrices, green - reduced friction

2 Soft Matrices, yellow - regular friction

2 Soft Matrices, red - high friction

REF 430 0516 0

Assortment

13 pieces

Vario-Soft 3 without integrated paralleling mandrel

2 vs 3 Patrices without paralleling mandrel

1 Matrix inserting instrument

2 Duplicating matrix

2 Wax matrix housing

2 Soft Matrices, green - reduced friction

2 Soft Matrices, yellow - regular friction

2 Soft Matrices, red - high friction

REF 430 0738 2

Assortment

13 pieces

Vario-Soft 3 Soft

2 vs 3 Patrices

1 Matrix inserting instrument

2 Duplicating matrix

2 Wax matrix housing

2 Soft Soft Matrices, green - reduced friction

2 Soft Soft Matrices, yellow - regular friction

2 Soft Soft Matrices, red - high friction

REF 430 0561 0