

Diamond grinding

- Diacryl Grinding Instruments dcs
- Diagen-Turbo-Grinder dtg
- Set-up grinding tool

- Special Diamonds for the Veneering Technique
- Diamond grinding tool dsl
- Diabolo / Diabolo Cleaner

- **FG-Diabolo**


FG-Diabolo

Fast, efficient processing of any type of zirconium oxide and hard dental materials.
Friction grip - first-class diamond grinders. Economic system of grinding tools with razor sharp diamonds, self-regenerating grit and extended durability. FG-Diabolo are sintered diamond grinding tools and are used in the turbine or in the T-hand-piece with FG adapter for grinding zirconium oxide and extremely hard materials. FG-Diabolo reduces processing times and renders working processes highly efficient.


FG adapter 1.6 to 2.35
for the use in the handpiece
Pack cont. 10 pieces
REF 340 0100 2
not included in the assortment

Assortment


6 pieces, 1 piece each
FG-Diabolo, fine grit
Bud, large
Flame
Cylinder, round head
Inverted cone
Torpedo
Bud, small
REF 330 0116 6

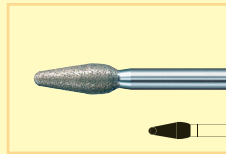


Bud, fine, large
REF FF 263 023

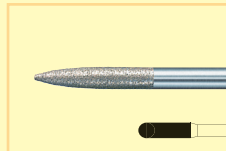


Flame, fine
REF FF 250 016

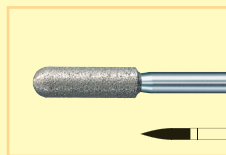


Cylinder, round head fine
REF FF 141 023

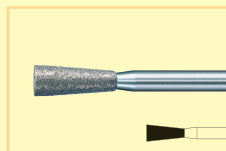

Inverted cone, fine
REF FF 227 023

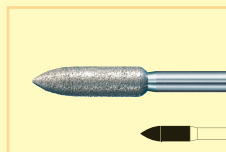

Torpedo, fine
REF FF 289 023

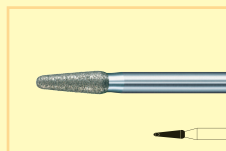

Bud, fine, small
REF FF 263 014



Accessories:



Diabolo Cleaner
1 piece
REF 340 0100 0

All the tools you need even when processing zirconium oxide

Indispensable tool for removing contaminations. Diabolo Cleaner guarantees constant cutting performance. Contaminated material is removed easily and quickly and new diamond crystals are exposed from the bronze binding material and integrated into the surface of the tool to enhance the cutting performance and reduce the grinding time.

Leading implant manufacturers recommend Diagen-Turbo-Grinders dtg. These abrasive tools have proved their suitability for reworking sintered zirconium frameworks and minimized the amount of work thanks to the special Diagen diamond binding material.



see also page 370

Sortiment

5 pieces, 1 piece each
Diagen-Turbo-Grinder dtg
REF 340 0020 0